

Non-Gynae Centrifuge Consumables

CellPath's non-gynaecological consumables provide solutions for both routine and large volume sample preparation. Our funnels are compatible with the EpreDia™ CytoSpin®*

**CellPath
Own Brand
Now Available
Worldwide**

Cytology Funnels and Clips

Manufactured to Our Exacting Standards

CellPath's Cytology funnels are the functional equivalent of the EpreDia™ CytoFunnel®* and are used to deposit a thin layer of cells in a clearly defined area of a microscope slide.

A pre-attached filter card comes with each funnel, which absorbs any excess fluid during the deposition of cells. Sealing caps are also included. Our range of Cytology funnels clips are compatible with the EpreDia™ CytoSpin®* centrifuge and come with easy-to-follow operational instructions.

Single Cytology Funnel with Filter Card – Disposable (0.4ml)

Accommodating sample volumes of up to 0.4ml, the single Cytology funnel provides lab staff with an easy, cost-efficient method of producing high-quality, thin-layer slide preparations.

The funnel comes with the option of having either a pre-attached white or brown filter card. The brown filter card absorbs samples at a slower rate than the white card, ideal for use with scanty specimens.

CellPath also offer individually wrapped single funnels.



KEY BENEFITS

- 1 Accommodates sample volumes of up to 0.4ml
- 2 Provides a consistent and homogeneous cell deposition area of 6mm
- 3 Filter cards come pre-attached
- 4 Optional brown filter card - absorbs samples at a slower rate than white
- 5 Brown filter card, excellent for scanty specimens such as Cerebral Spinal Fluid
- 6 Funnel sealing caps included

Item Code	Description	Quantity
NND-4000-09A	Disposable Cytology Funnel Single 0.4ml - Individually Wrapped	50
NND-4000-10A	Disposable Cytology Funnel Single 0.4ml - Bulk	500
NND-4000-10B	Disposable Cytology Funnel Single 0.4ml - Individually Wrapped	500
NND-4001-09A	Disposable Cytology Funnel (Brown) Single 0.4ml - Individually Wrapped	50
NND-4001-10A	Disposable Cytology Funnel (Brown) Single 0.4ml - Individually Wrapped	500

Double Cytology Funnel with Filter Card – Disposable (0.25ml)

The double Cytology funnel can hold sample volumes of up to 0.25ml and is designed to process twice the number of sample preparations per run.

As a result, it produces two 6mm cell deposition areas on each slide. Like the single funnel, they are also available individually wrapped.



KEY BENEFITS

- 1 Accommodates sample volumes of up to 0.25ml
- 2 Provides two consistent and homogeneous cell deposition areas of 6mm per slide
- 3 Funnel markings correspond to slide markings for simple identification
- 4 Pre-attached filter card
- 5 Funnel sealing caps included


Item Code	Description	Quantity
NND-4004-08A	Disposable Cytology Funnel Double 0.25ml - Individually Wrapped	25
NND-4004-10A	Disposable Cytology Funnel Double 0.25ml - Bulk	500
NND-4004-12A	Disposable Cytology Funnel Double 0.25ml - Individually Wrapped	250

Cytology Funnel Clips

CellPath provide a stainless steel clip, which holds the Cytology funnel snugly against the microscope slide during the centrifugation process.

The stainless steel clips can be sanitised after each use.

The clip is compatible with the single and double chamber Cytology funnels.

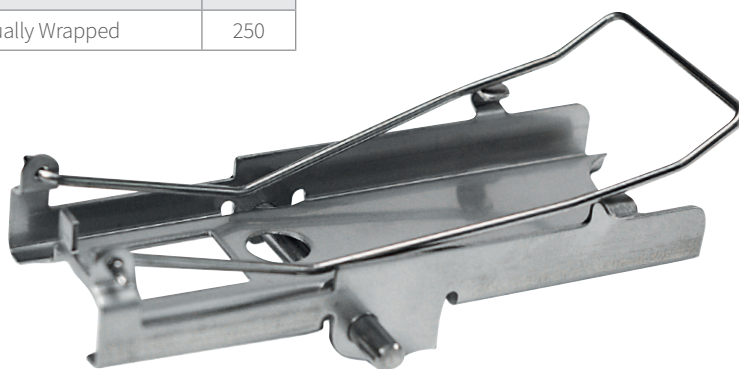

KEY BENEFITS

- 1 Compatible with both CellPath's single and double chamber Cytology funnels

Item Code	Description	Quantity
NND-4008-50A	Stainless Steel Cytology Funnel Clip	6

The Cytology Funnels and Clips are fully compatible with our range of Cytology Slides.



CellPath Ltd

80 Mochdre Enterprise Park,
Newtown, Powys, SY16 4LE, UK

T: +44 (0)1686 611333

E: sales@cellpath.com

cellpath.com



*CytoSpin® and CytoFunnel® are a registered trademark of EpreDia™